

SCAN 50 Manual colony counter

# A galaxy of colonies with a clear view

SCAN 50 interscience




**Scan® 50** is a new generation colony counter, offering contrasted colony lighting while protecting the user's eyes. It improves the ergonomics of colony counting and prevents MusculoSkeletal disorders (MSD). Traceability is ensured by USB port.



# **Contrasting colony lighting:** Dark Field 2

Dark Field 2 is an indirect lighting system with adjustable intensity. Light reflectors diffuse the light at all points without any reflections or shadows. It makes all colonies stand out on the surface, in the mass and around the edge of the dish.

Flashes of light are eliminated and the **operator's eyes are protected.** 

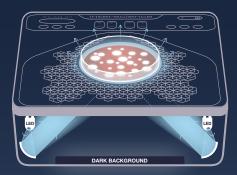





## **Easy disinfection**

To keep the workstation clean, the glass and the smooth handrest are removable.

Easy to disinfect, they ensure **perfect** hygiene of the workstation.

Headphones jack



Scan<sup>®</sup> 50 pro



SCAN 50 PTO

interscience

## **Ergonomics redesigned**

#### **Audio Comfort**

• 9 different sounds and headphones jack to work without disturbing colleagues


HandPad

#### Light signal

• Light signal for each counted colony

#### **Optimal working posture**

- Inclined reading surface with magnifying lens<sup>(2)</sup>, hand rest and HandPad(1)
- Prevents Muscular-Skeletal Disorders

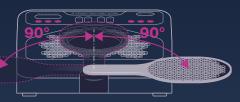


#### Ambidextrous

- Duplicated count buttons on both sides
- HandPad, 180° arm-rest that turns on the machine

## **Technical specifications** & Accessories

Delivered with: 1 power cord, 1 USB cable, 1 black marker pen, 1 white diffusing background, 1 user manual, 1 headphones, (only on Scan® 50 pro), 1 conformity certificate / warranty card. The magnifying lens is optional (Ref. 435 001)

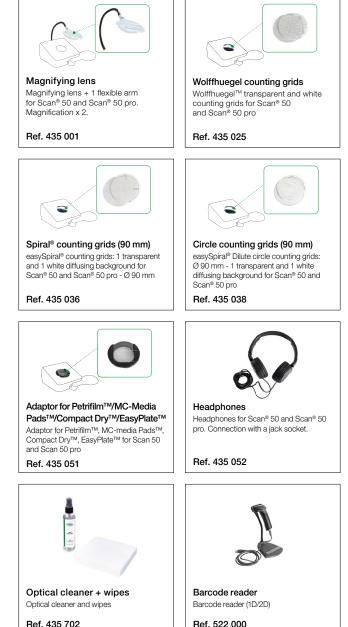

Ref. 435 702





		Scan <sup>®</sup> 50	Scan <sup>®</sup> 50 pro
OVERVIEW	Reference	435 050	435 055
	Painted stainless steel shell	✓	✓
	LIMS/SIL connection	✓	✓
	USB connection	✓	✓
	Counting on pour, surface, Spiral® and circle plated Petri dishes	✓	✓
	Counting on chromogenic dishes	✓	✓
	Counting on Petrifilm™, Compact Dry™, MC-Media Pads™, EasyPlate™, filtration membranes	✓	✓
	White LED lighting without reflections	✓	✓
	HandPad: ergonomic arm-rest	-	✓
	Headphones included	-	✓
COUNTING	Counting	Manual (0-1999 CFU)	
	Manual control to add or substract colonies	✓	✓
	Counting light signal	✓	✓
	Adjustable counting beeps	9 different sounds / 5 volume levels	
	Counting zone	Touch surface with adjustable sensitivity	
SPECIFICATIONS	White LED Lighting technology	Dark Field 2	
	LED Lighting system	5 intensity levels	
	Petri dish size	From Ø 55 to Ø 90 mm	
	USB Data export	To LIMS, as CSV spreadsheet (Excel™) or text file (Word™)	
	Voltage - Frequencies	100-240 V~ 50-60 Hz	
	Warranty	3-year (after recording the warranty card)	
	In compliance with	ISO 7218 and FDA BAM (Bacteriological Analytical Manual)	
WEIGHT AND DIMENSIONS	Dimensions (w x d x h): Scan	27.3 x 29.5 x 15.1 cm	37.6 x 35.2 x 15.1 cm (HandPad open)
	Net weight	4.05 kg	4.45 kg
	Box (w x d x h)	35 x 45 x 40 cm	35 x 45 x 40 cm
	Gross weight: Scan	5.50 kg	6.15 kg
WEEE 2002/96/EC ROHS 2002/95/EC FOR CE UK Products made for INTERSCIENCE by Interlab, an ISO 9001 certified company.			

Contact us for IQOQ and warranty extension

#### interscience

#### PARIS

Phone: +33 (0)1 34 62 62 61 - Email: info@interscience.com

#### FRANKFURT

Phone: +49 611 7238 7770 - Email: sales.germany@interscience.com

#### BOSTON

Phone: +1 781 937 0007 - Email: sales.usa@intersciencelab.com

SHANGHAI Phone: +86 (0)21-64739390 - Email: sales.china@interscience.cn

#### SINGAPORE

Phone: +65 6977 7232 - Email: sales.asia@interscience.com TOKYO

Phone: +81 3 6712 9715 - Email: sales.japan@interscience.com